everStick®C&B

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951327-NU7708
30000523

DISTRIBUTED BY
GC CORPORATION
76-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8585, Japan

GC AMERICA INC.
3737 West 127th Street, Alsip, IL 60803 U.S.A.
TEL: +1-708-597-0900

GC ASIA DENTAL PTE. LTD.
5 Tampines Central 1, #06-01 Tampines Plaza, Singapore 529541
TEL: +65 6546 7588

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Temperature limit
2 °C / 35 °F

8 °C / 46 °F

Caution: US Federal law restricts this device to sale by or on the order of a dentist.


everStick®C&B

EN FIBRE REINFORCEMENT FOR COMPOSITE BRIDGES
BG ФИБРОВЛАКНА ЗА УКРЕПВАНЕ НА КОМПОЗИТНИ МОСТОВЕ
CS VLÁKNOVÉ VÝZTUŽE PRO KOMPOZITNÍ MŮSTKY
HR VLAKNA ZA POJAČANJE ZA KOMPOZITNE MOSTOVE
HU ŰVEGSZÁLAS MEGERŐSÍTÉS KOMPOZIT HIDAKHOZ
PL WZMOCNIENIE Z WŁÓKNA SZKLANEGO DO MOSTÓW KOMPOZYTOWYCH
RO FIBRE CONSOLIDANTE PENTRU PUNȚI DIN COMPOZIT
RU АРМИРУЮЩЕЕ СТЕКЛОВОЛОКНО ДЛЯ МОСТОВИДНЫХ КОНСТРУКЦИЙ ИЗ КОМПОЗИТОВ
SK VÝSTUŽ ZO SKLENÝCH VLÁKIEN PRE VLÁKNAMI VYSTUŽENÉ KOMPOZITNÉ MOSTÍKY
SL VLAKNATE OJAČITVE ZA KOMPOZITNE MOSTIČKE
SR FIBER OJAČANJE ZA KOMPOZITNE MOSTOVE
UK ПІДСИЛЕНИЙ СКЛОВОЛОКНОМ МАТЕРІАЛ ДЛЯ КОМПОЗИТНИХ МОСТІВ
TR KOMPOZIT KÖPRÜLER İÇİN FİBER GÜÇLENDİRME
LV ŠKIEDRU STIPRINĀJUMS KOMPOZĪTA TILTIEM
LT KOMPOZITINIUS TILTUS SUTVIRTINANČIOS SKAIDULOS
ET FIBERTUGEVDUS KOMPOSIITSILDADELE

GC

everStick®C&B 1. A surface retained bridge in the anterior region

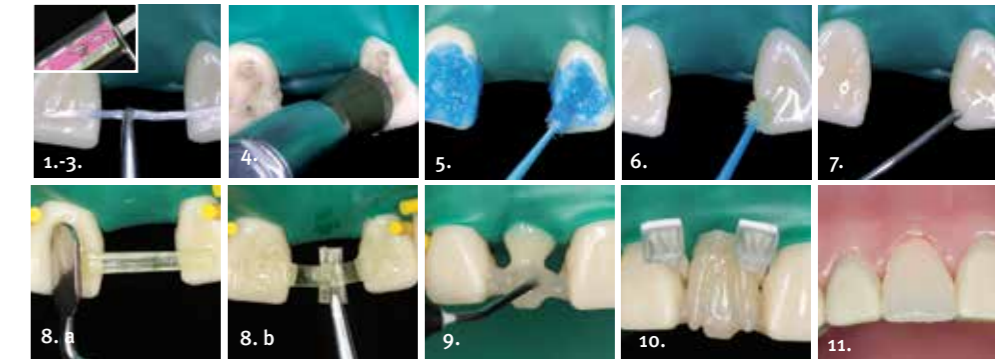EN	FIBRE REINFORCEMENT FOR COMPOSITE BRIDGES	7	SK	VÝSTUŽ ZO SKLENÝCH VLÁKIEN PRE VLÁKNAMI VYSTUŽENÉ KOMPOZITNÉ MOSTÍKY	97
BG	ФИБРОВЛАКНА ЗА УКРЕПВАНЕ НА КОМПОЗИТНИ МОСТОВЕ	18	SL	VLAKNATE OJAČITVE ZA KOMPOZITNE MOSTIČKE	107
CS	VLÁKNOVÉ VÝZTUŽE PRO KOMPOZITNÍ MŮSTKY	30	SR	FIBER OJAČANJE ZA KOMPOZITNE MOSTOVE	116
HR	VLAKNA ZA POJAČANJE ZA KOMPOZITNE MOSTOVE	40	UK	ПІДСИЛЕНИЙ СКЛОВОЛОКНОМ МАТЕРІАЛ ДЛЯ КОМПОЗИТНИХ МОСТІВ	126
HU	ÜVEGSZÁLAS MEGERŐSÍTÉS KOMPOZIT HIDAKHOZ	50	TR	KOMPOZIT KÖPRÜLER İÇİN FIBER GÜÇLENDİRME	138
PL	WZMOCNIENIE Z WŁÓKNA SZKLANEGO DO MOSTÓW KOMPOZYTOWYCH	60	LV	ŠĶĪEDRU STIPRINĀJUMS KOMPOZĪTA TILTIEM	148
RO	FIBRE CONSOLIDANTE PENTRU PUNȚI DIN COMPOZIT	72	LT	KOMPOZITINIUS TILTUS SUTVIRTINANČIOS SKAIDULOS	157
RU	АРМИРУЮЩЕЕ СТЕКЛОВОЛОКНО ДЛЯ МОСТОВИДНЫХ КОНСТРУКЦИЙ ИЗ КОМПОЗИТОВ	83	ET	FIIBERTUGEVDUS KOMPOSIITSILDADELE	166

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- 1 pontic (2-unit bridge): 2 everStickC&B fibre bundles

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- 1 pontic (2-unit bridge): 3 everStickC&B fibre bundles
- Maximum amount of pontics is 1.

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1. Use articulating paper to **ensure that there is enough space in the occlusion**

- for making a surface-retained bridge.
- Measure** the length of the fibre required to prepare an everStickC&B fibre frame using, for example, a periodontal probe or dental floss. The fibre should cover approximately two-thirds of the width of the supporting tooth's crown. Open the foil package and use tweezers to pull only the required amount of silicone bedding out of the package.
Note: Use powder-free gloves when handling everStickC&B fibre.
 - Cut** the required amount of fibre together with its silicone bedding using sharp scissors. Shield the piece of fibre from light by placing it under a cover during preparation of the teeth to be bonded. Close the foil bag tightly with its sticker. Store the bag in a refrigerator (2-8°C, 35-46°F) when you do not need it.
 - Clean** the area to be bonded using a pumice and water mix, rinse with water, and air-dry the area.
 - Etch** the surfaces of the teeth in the bonding areas with orthophosphoric acid according to the bonding agent manufacturer's instructions. The recommended enamel etching time for surface-retained areas is 45 to 60 seconds. The area to be etched must be sufficiently large. Preferably etch a slightly wider area than needed rather than too small. Careful etching and bonding techniques ensure a reliable bonding of the bridge to the teeth. Rinse with water and air-dry the tooth surfaces carefully after etching.
Note: Keep the working area dry during bonding and layering of composite. Use of a rubber dam is highly recommended to keep working conditions dry.
 - Bond** using the composite bonding technique described in the manufacturer's instructions of your bonding agent. Apply the bonding agent to the entire area to be bonded. Light-cure the bonding agent as recommended by the manufacturer.
 - Apply a thin layer of flowable composite** (for example G-ænial Universal Flo) to the bonded surfaces of the teeth. **Do not light-cure the flowable composite during this phase.**
 - Place the fiber.** Remove the white cover paper and use tweezers to pick the fibre up from the silicone groove. Place the fibre bundle on the tooth on top of the uncured, flowable composite. Press the other end of the fibre bundle or the entire fibre bundle tightly onto the surface of the tooth using a StickREFIX D silicone instrument or a StickSTEPPER hand instrument (in case of

cavity-retained bridge, use a Stick CARRIER hand instrument). Light-cure for 5 to 10 seconds, but at the same time, protect the other end of the fibre bundle from premature light-curing using a wide StickSTEPPER instrument. Position and light-cure the rest of the fibre one tooth/pontic area at a time, as above. Only light-cure the fibre for 5 to 10 seconds per tooth at this stage.

Note:

- a) Spread the fibre bundle wider on the surface-retained areas to create a more extensive bonding area.
- b) Do not place the fibre too close to the gingiva approxinally, so that the cleaning spaces are not covered.
- c) Attach the fibre frame as incisally as possible. This allows maximum support for the bridge in the anterior region.

- d) The fibre in the pontic section is meant to run along the middle of the future pontic in the lingual-labial direction.
- e) When necessary, add transverse fibre to support the composite pontic. The transverse fibre should be placed on the occlusal side of the fibre frame.

9. **Cover & Light-cure the fiber frame.**

After the initial light-curing, cover the entire fibre frame with a thin layer of composite. Light-cure the entire fibre frame for 40 seconds, one unit at a time.

Note: The fibres must be covered entirely with composite, including the interproximal areas. However, there should be enough space to allow the patient to clean the bridge and the approximal areas.

10. **Layer the pontic** using composite according to the composite manufacturer's instructions. If you do not use a rubber dam,

you can use, for example, a rubber dam strip, a sectional matrix, or a plastic strip as a moisture barrier against the gingiva. Build the core and base using flowable composite. This simplifies formation of the base of the pontic. The pontic region next to the gingiva should have a light contact point, and the form should be self-cleaning. Build the dentine parts of the tooth with dentine shades and the enamel parts with enamel shades. If you wish, you can also use characterisation shades.

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Note:

Use powder-free gloves when handling everStickC&B fibre.

2. **Cut** the required amount of fibre together with its silicone bedding, using sharp scissors. Shield the piece of fibre from light by placing it under a cover during preparation of the teeth to be bonded. Close the foil bag tightly with its sticker. Store the bag in a refrigerator (2-8°C, 35-46°F) when you do not need it.
3. **Clean** the unprepared areas of the teeth to be bonded using pumice and water, rinse with water, and air-dry the area.
4. **Etch** the surfaces of the teeth in the bonding

areas with ortho-phosphoric acid according to the bonding agent manufacturer's instructions. The recommended enamel etching time for surface-retained areas is 45 to 60 seconds. The area to be etched must be sufficiently large.

Preferably etch a slightly wider area than needed rather than too small. Careful etching and bonding techniques ensure a reliable bonding of the bridge to the teeth. Rinse with water and air-dry the tooth surfaces carefully after etching.

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5. **Bond** using the composite bonding technique described in the manufacturer's instructions of the bonding agent. Apply the bonding agent to the entire area to

be bonded. Light-cure the bonding agent as recommended by the manufacturer.

6. **Apply a thin layer of flowable composite** (for example, G-aenial Universal Flo) to the bonded surfaces of the teeth and/or the floor of the cavity. Do not light-cure the composite during this phase.
7. **Place the fibre.** Remove the white cover paper and use tweezers to pick the fibre up from the silicone groove. Place the fiber bundle on the tooth on top of the uncured, flowable composite. Press the other end of the fibre bundle or the entire fibre bundle tightly onto the surface of the tooth using a StickREFIX D silicone instrument, a StickSTEPPER hand instrument, or into the cavity using a StickCARRIER hand instrument. Light-cure for 5 to 10 seconds, but at the same time protect the other end of the fibre

bundle from pre-mature light-curing using a wide StickSTEPPER instrument. Position and light-cure the rest of the fibre one tooth/pontic at a time, as above. Only light-cure the fibre for 5 to 10 seconds per tooth at this stage.

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- e) Position and light-cure any additional fibres as described above. Use flowable composite to attach the fibres together.
 - f) When necessary, add transverse fibres to support the composite pontic.
8. **Cover and Light-Cure the fibers.** After positioning and light-curing, cover the entire fibre frame with a thin layer of composite. Light-cure the entire fibre frame for 40 seconds, one unit at a time.
- Note:** The fibres must be covered entirely by composite, including the interproximal areas. However, there should be enough space to allow the patient to clean the bridge and the approximal areas.
9. **Layer the pontic** using composite according to the composite manufacturer's instructions. If you do not use a rubber

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2. Remove any temporary restorations and check that the work fits.
3. Prosthetic pre-treatment
 - a) Use a carborundum bur to lightly roughen the surfaces to be cemented. Rinse with water and air-dry the surfaces. **Note:** Do not use sandblasting with StickNET fibre.
 - b) Apply enamel bonding agent (for example, StickRESIN) to the roughened bonding surfaces** to activate it, protect it from light and let it take effect for 3 to 5 minutes (you can use, for instance, a metal cup as protection against light). Carefully remove excess bonding agent

with air, because an overly thick layer of bonding agent prevents the work from fitting perfectly. Light-cure the bonding agent for 10 seconds before cementation.

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Preparing the teeth

4. Clean the surface-retained areas using a pumice and water mix.
5. Etch the surfaces of the teeth over a wide area according to the cement manufacturer's instructions. The recommended enamel etching time for surface-retained areas is 45 to 60 seconds using a 37%

orthophosphoric acid. Rinse with water and air-dry the tooth surfaces thoroughly.

- Bond the teeth according to the cement manufacturer's instructions.

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Cementing

- Apply a dual cure or chemical cure composite cement onto the cementing surfaces of the work, and position the work in place.

Note: Use dual cure or chemical cure composite cements for cementing fibre work. Phosphate and glass ionomer cements are NOT suitable for cementing fibre work.

- Remove excess cement and apply oxygen-blocking gel (for instance, glycerol gel) onto the marginal areas.
- Light-cure dual cure cement according to

the cement manufacturer's instructions.

- Check and adjust the occlusion. Finish. Be careful not to cut the fibres when finishing the approximal areas.

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- Before preparation or extraction of the tooth/ teeth, take an impression of the working area.
- Prepare the abutment teeth. To prepare the everStickC&B fibre frame, measure the length of fibre needed with, for example, dental floss.
- Cut the required length of fibre together

with its silicone bedding, using sharp scissors. Shield the piece of fibre from light by placing it under a cover until it will be used. Close the foil bag tightly with its sticker, and store the bag in a refrigerator (2-8°C, 35- 46°F).

4a-4d

Apply a thin layer of light curing adhesive* (e.g. StickRESIN) on top of the abutments. Do not light-cure the adhesive yet. Place the fibre bundle on top of the abutments. Use a StickREFIX D silicone instrument to press the fibre into the correct position and shape.

To achieve maximum reinforcement, bend the fibre in the pontic area to a curved form as close as possible to the gingiva. However, leave approximately 1-2 mm of space between the fibre and the gingiva for the temporary crown and bridge

material. Aproximally do not place the fibre too close to the gingiva, so as not to cover the cleaning spaces. Initially light-cure the fibre through the silicone for 10 seconds. Finalise the polymerisation by light-curing the entire fibre frame for 40 seconds per area.

*** NOTE:** The light curing adhesive used must be monomer-based and must not contain etchants or solvents (acetone, alcohol, water, etc.). Do not etch the abutments or use dentin primer prior to adhesive application.

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5. Follow the user instructions of the temporary crown and bridge material manufacturer. Dispense temporary material first between the bridge frame and the gingiva, to eliminate air gaps. Fill the impression

with temporary material, place it over the prepared teeth, and allow it to set.

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STORING: everStick products should always be stored in a refrigerator (2-8°C, 35-46°F). In addition, the products should be protected from light by packing them in the sealed foil package after use. An elevated temperature and exposure to bright light may shorten the lifetime of everStick products.

Prior to application, the products are taken

out of the refrigerator and the foil package opened, but kept away from bright daylight or artificial light. While cutting the fibre bundle, the rest of the fibre bundle inside the foil package should be kept covered from light. Immediately after cutting a sufficient length for the fibre construction, the foil package is carefully resealed and returned to the refrigerator.

(Shelf life: 2 years from date of manufacture)

PACKAGES

everStick Intro

8 cm everStickC&B; 8 cm everStickPERIO;
1x StickREFIX D silicone instrument

everStick COMBI

8cm everStickC&B; 8cm everStickPERIO;
30 cm² everStickNET; 5x everStickPOST Ø 1.2;

1x StickSTEPPER hand instrument;
1x StickREFIX D silicone instrument"

everStick Starter Kit:

8 cm everStickC&B; 8 cm everStickPERIO;
30 cm² everStickNET; 5 x everStickPOST
Ø 1.2; 5 ml StickRESIN bottle; 2 ml G-æniel
Universal Flo syringe; 20 dispensing tips,
1 light protective cover; 1x StickREFIX D
silicone instrument; 1x StickSTEPPER;
1x StickCARRIER

Refills

2 x 12 cm fibre bundle
1 x 8 cm fibre bundle

NOTE: The everStick products should be used clinically with care and the patient should be warned not to abrade the fitting surface so as to avoid exposing irritation-causing fibres.

If the surface of the fibre bundle feels dry, but it is fully bendable and not polymerized, adding a drop of resin (such as StickRESIN) will return the flexibility/workability of the material.

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WARNING: Unpolymerised resin can cause skin sensitisation to acrylates in some people. If your skin comes in contact with resin, wash it thoroughly with soap and water. Avoid contact of uncured material with skin, mucous membrane, or eyes. Unpolymerised everStick products may have a slight irritating effect and lead to sensitization to methacrylates in rare cases. The use of powder free gloves is recommended with everStick products. Polymerize everStick before waste disposal.

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Viimati parandatud: 07/2019

EN Undesired effects - Reporting:

If you become aware of any kind of undesired effect, reaction or similar events experienced by use of this product, including those not listed in this instruction for use, please report them directly through the relevant vigilance system, by selecting the proper authority of your country accessible through the following link: https://ec.europa.eu/growth/sectors/medical-devices/contacts_en as well as to our internal vigilance system: vigilance@gc.dental

In this way you will contribute to improve the safety of this product. **BG** Докладване на нежелани ефекти:

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https://ec.europa.eu/growth/sectors/medical-devices/contacts_en

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G-ænia!®Universal Flo

everStick®C&B

Fibre reinforcement for composite bridges

DISTRIBUTED BY

GC CORPORATION
76-1 Hasunuma-cho, Itabashi-ku,
Tokyo 174-8585, Japan

GC EUROPE N.V.
Researchpark Haasrode-Leuven 1240,
Interleuvenlaan 33,
B-3001 Leuven, Belgium
TEL: +32 16 74 10 00

GC AMERICA INC.
3737 West 127th Street, Alsip, IL 60803 U.S.A.
TEL: +1-708-597-0900
www.gcamerica.com

RESPONSIBLE MANUFACTURER IN CANADA

GC AMERICA INC.
3737 W. 127th Street, Alsip, IL 60803 U.S.A.

GC ASIA DENTAL PTE. LTD.
5 Tampines Central 1
#06-01 Tampines Plaza
Singapore 529541 TEL: +65 6546 7588

GC AUSTRALASIA DENTAL PTY. LTD.
1753 Botany Rd, Banksmeadow,
NSW 2019, Australia
TEL: +61 2 9301 8200

GC SOUTH AMERICA
Rua Heliadora, 399,
Santana - São Paulo, SP, Brasil
CEP: 02022-051
TEL: +55-11-2925-0965
CNPJ: 08.279.999/0001-61
RESP. TÊC: MAYARA DE SANTIS RIBEIRO -
CRO/SP 105.982

Fibre type: Silanated E-Glass fibre impregnated with Bis-GMA and PMMA

Form: Unidirectional fibre bundle
Diameter: ~1.5-1.7 mm

Keep away from sunlight

Temperature limit

Caution: US Federal law restricts this device to sale by or on the order of a dentist.
NU7610 - 30000524

951328 – NU7610
30000524


everStick®C&B

EN
DE
FR
IT
ES
NL
DA
SV
PT
EL
FI
ZH
ID






FIBRE REINFORCEMENT FOR COMPOSITE BRIDGES
GLASFASERVERSTÄRKUNG FÜR COMPOSITEBRÜCKEN
FIBRE DE RENFORT POUR BRIDGES COMPOSITE
RINFORZO IN FIBRA PER PONTI IN COMPOSITO
FIBRA DE REFUERZO PARA PUENTES DE COMPOSITE
VEZELVERSTERKING VOOR COMPOSITBRUGGEN
FIBERFORSTÆRKNING TIL KOMPOSITBROER
FIBERFÖRSTÄKNING TILL KOMPOSITBROAR
FIBRA PARA REFORÇO DE PONTES DE RESINA COMPOSTA
ΙΝΕΣ ΥΑΛΟΝΗΜΑΤΩΝ ΓΙΑ ΕΝΙΣΧΥΣΗ ΓΕΦΥΡΩΝ ΑΠΟ ΠΟΛΥΜΕΡΗ ΥΛΙΚΑ
LASIKUITULUJITE KUITULUJITTEISIIN YHDISTELMAMUOVISILTOIHIN
用于各种合成牙桥制作的复合树脂纤维
SERAT YANG DIPERKUAT UNTUK JEMBATAN KOMPOSIT

everStick®
C&B 1. A surface retained bridge in the anterior region

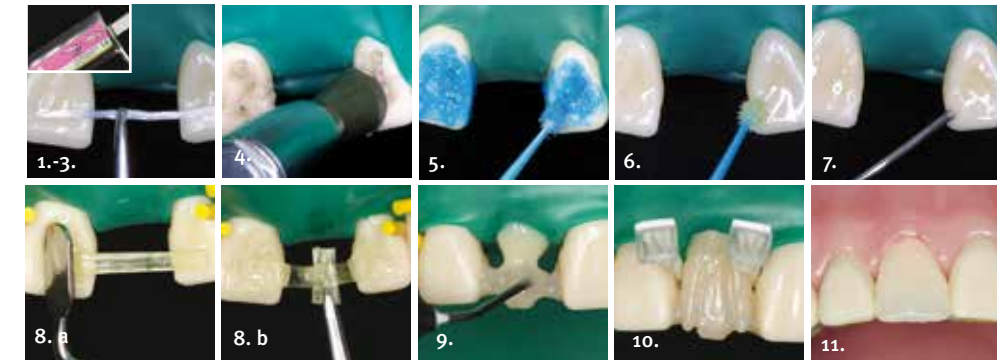everStick[®] C&B 4. Direct fibre-reinforced temporary bridge with everStickC&B

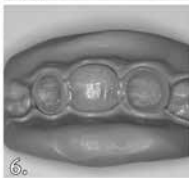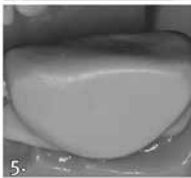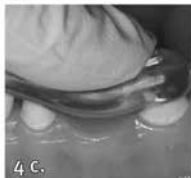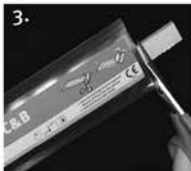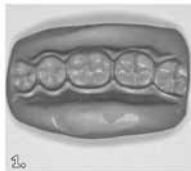EN	Fibre reinforcement for composite bridges	7	SV	Fiberförstärkning till kompositbroar	93
DE	Glasfaserverstärkung für Compositebrücken	18	PT	Fibra para reforço de pontes de resina composta	104
FR	Fibre de renfort pour bridges composite	31	EL	Ίνες υαλονημάτων για ενίσχυση γεφυρών από πολυμερή υλικά	116
IT	Rinforzo in fibra per ponti in composito	43	FI	Lasikuitulujite kuitulujitteisiin yhdistelmämuovisiltoihin	130
ES	Fibra de refuerzo para puentes de composite	56	ZH	用于各种合成牙桥制作的复合树脂纤维	141
NL	Vezelversterking voor composietbruggen	69	ID	Serat yang diperkuat untuk jembatan komposit	149
DA	Fiberforstaerkning til kompositbroer	82			

WHAT IS everStickC&B?

everStickC&B is a fibre reinforcement made from glass fibres that add strength to composites.

everStickC&B fibre reinforcement is made from unidirectional glass fibres and a polymer/resin gel matrix. The polymer/resin gel holds the individual glass fibres in a bundle, which facilitates handling of the fibre bundle. The fibre bundle is flexible and tacky, which allows it easily and reliably bond to teeth.

INDICATIONS

everStickC&B fibre reinforcement is recommended for use in the following applications:

- Inlay, onlay, full coverage and hybrid bridges
- Surface-retained bridges
- Temporary bridges
- Bridges made indirectly at the chairside or

in a dental laboratory

- Reinforcement and repair of removable dentures and appliances

CONTRA-INDICATION

In rare cases the product may cause sensitivity in some people. If any such reactions are experienced, discontinue the use of the product and refer to a physician.

COMPATIBLE MATERIALS FOR everStickC&B REINFORCEMENT

- methacrylate-based dental composites (light, chemical and dual cure)
- methacrylate and acrylic resins/monomers, and polymerizable bonding adhesives
- methacrylate-based dental composite resin cements (light, chemical and dual cure)
- denture base acrylics

COMPATIBLE MATERIALS FOR REPAIRING everStickC&B RESTORATIONS

- solvent free adhesive resins
- methacrylate-based resins
- denture base acrylics

INSTRUCTIONS FOR USE

Frame Design

The everStickC&B fibre frame can be surface-retained on the bridge's abutment teeth or by embedding the glass fibre in prepared cavities. A combination structure, containing both surface retention and cavity preparation, provides the best support in structures supporting themselves via molar, premolar or canine teeth. The fibre frame can be attached buccally, lingually and/or occlusally, depending on the clinical situation. The fibre frame should also be bonded to the teeth for their entire length, including interproximal areas.

Prior to fibre placement, articulating paper should be used to check whether there is enough room for the fibre frame in occlusion, particularly in surface-retained constructions. At the site of occlusal contact, the optimal height of veneering composite to be layered on top of the fibre frame is 1-2 mm so that it will not detach from the fibre. This must be taken into consideration when planning the fibre frame design and the required preparations.

Placement and number of fibres

The fibres must always be placed at the location where there is most tensile stress in the structure.

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The fibre should be placed as incisally as possible and should run in the middle of the pontic towards the labial surface.

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The fibre should be placed opposite to the occlusal surface of the pontic, close to the gingiva. In the pontic area, a wide fibre frame gives more support to the composite pontic. For the best support, especially in large pontic constructions, add short transverse fibres to the fibre frame to support the pontic.

Number of pontics and fibres in everStickC&B fixed bridges:

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- 1 pontic (3-unit bridge): 1 everStickC&B fibre bundle
- 2 pontics (4-unit bridge): 2 everStickC&B fibre bundles
- 3 pontics (5-unit bridge): 3 everStickC&B fibre bundles

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- 1 pontic (3-unit bridge): 2 everStickC&B fibre bundles
- 2 pontics (4-unit bridge): 3 everStickC&B fibre bundles
- 3 pontics (5-unit bridge): 4 everStickC&B fibre bundles

Maximum amount of pontics is 3.

Number of pontics and fibres in everStickC&B cantilever bridges:

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- 1 pontic (2-unit bridge): 2 everStickC&B fibre bundles

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- 1 pontic (2-unit bridge): 3 everStickC&B fibre bundles

Maximum amount of pontics is 1.

INSTRUCTIONS BY INDICATION

I. SURFACE-RETAINED BRIDGE WITH everStick C&B FIBRE IN THE ANTERIOR REGION

A bridge in the anterior region is made using one everStickC&B fibre per pontic. The use of a rubber dam is highly recommended to keep the working area dry.

1. Use articulating paper to **ensure that there is enough space in the occlusion** for making a surface-retained bridge.

- Measure** the length of the fibre required to prepare an everStickC&B fibre frame using, for example, a periodontal probe or dental floss. The fibre should cover approximately two-thirds of the width of the supporting tooth's crown. Open the foil package and use tweezers to pull only the required amount of silicone bedding out of the package.
Note: Use powder-free gloves when handling everStickC&B fibre.
- Cut** the required amount of fibre together with its silicone bedding using sharp scissors. Shield the piece of fibre from light by placing it under a cover during preparation of the teeth to be bonded. Close the foil bag tightly with its sticker. Store the bag in a refrigerator (2-8°C, 35,6-46,4°F) when you do not need it.
- Clean** the area to be bonded using a pumice and water mix, rinse with water, and air-dry the area.
- Etch** the surfaces of the teeth in the bonding areas with orthophosphoric acid according to the bonding agent manufacturer's instructions. The recommended enamel etching time for surface-retained areas is 45 to 60 seconds. The area to be etched must be sufficiently large. Preferably etch a slightly wider area than needed rather than too small. Careful etching and bonding techniques ensure a reliable bonding of the bridge to the teeth. Rinse with water and air-dry the tooth surfaces carefully after etching.
Note: Keep the working area dry during bonding and layering of composite. Use of a rubber dam is highly recommended to keep working conditions dry.
- Bond** using the composite bonding technique described in the manufacturer's instructions of your bonding agent. Apply the bonding agent to the entire area to be bonded. Light-cure the bonding agent as recommended by the manufacturer.
- Apply a thin layer of flowable composite** (for example G-aenial Universal Flo) to the bonded surfaces of the teeth. **Do not light-cure the flowable composite during this phase.**
- Place the fiber.** Remove the white cover paper and use tweezers to pick the fibre up from the silicone groove. Place the fibre bundle on the tooth on top of the uncured, flowable composite. Press the other end of the fibre bundle or the entire fibre bundle tightly onto the surface of the tooth using a StickREFIX D silicone instrument or a StickSTEPPER hand instrument (in case of cavity-retained bridge, use a Stick CARRIER

hand instrument). Light-cure for 5 to 10 seconds, but at the same time, protect the other end of the fibre bundle from premature light-curing using a wide StickSTEPPER instrument. Position and light-cure the rest of the fibre one tooth/pontic area at a time, as above. Only light-cure the fibre for 5 to 10 seconds per tooth at this stage.

Note:

- a) Spread the fibre bundle wider on the surface-retained areas to create a more extensive bonding area.
- b) Do not place the fibre too close to the gingiva approxinally, so that the cleaning spaces are not covered.
- c) Attach the fibre frame as incisally as possible. This allows maximum support for the bridge in the anterior region.
- d) The fibre in the pontic section is meant to

run along the middle of the future pontic in the lingual-labial direction.

- e) When necessary, add transverse fibre to support the composite pontic. The transverse fibre should be placed on the occlusal side of the fibre frame.
9. **Cover & Light-cure the fiber frame.** After the initial light-curing, cover the entire fibre frame with a thin layer of composite. Light-cure the entire fibre frame for 40 seconds, one unit at a time.
- Note:** The fibres must be covered entirely with composite, including the interproximal areas. However, there should be enough space to allow the patient to clean the bridge and the approximal areas.
10. **Layer the pontic** using composite according to the composite manufacturer's instructions. If you do not use a rubber dam, you can use, for example, a

rubber dam strip, a sectional matrix, or a plastic strip as a moisture barrier against the gingiva. Build the core and base using flowable composite. This simplifies formation of the base of the pontic. The pontic region next to the gingiva should have a light contact point, and the form should be self-cleaning. Build the dentine parts of the tooth with dentine shades and the enamel parts with enamel shades. If you wish, you can also use characterisation shades.

Note:

- a) Under the occlusal contact, the optimal thickness of the composite to be layered on top of the fibre frame is 1-2 mm so that the composite will not fracture from the fibre.
 - b) Remember to preserve cleaning spaces.
11. **Finish the bridge and adjust it into occlusion.**

Note:

- a) If, after placing the fibre, you notice that the fibre is too long, shorten it using a diamond bur during the finishing phase of the bridge. Apply enamel bonding agent (for example, StickRESIN) to the exposed fibre surface to activate it and carefully remove excess bonding agent with air. Light-cure the bonding agent for 10 seconds and carefully cover the fibre again with composite.
- b) Be careful not to cut the fibres in the finishing phase

II BRIDGE USING everStickC&B FIBRE IN THE POSTERIOR REGION

Frame design

In order to prepare a bridge in the posterior region, you need two everStickC&B fibres for bridges with one pontic. Use three

everStickC&B fibres for bridges with two pontics. Add one or two short, transverse fibres to the fibre frame to support the composite in the pontic. The use of rubber dam is highly recommended to ensure optimal working conditions. The everStickC&B fibre frame can be surface-retained on the bridge's abutment teeth and/or the glass fibre reinforcement can be embedded within prepared cavities. A combination structure, containing both a surface-retained and an inlay retainer, gives the best strength in structures supporting themselves on molar, premolar or canine teeth. By removing old restorations, you can obtain space for the fibre frame without additional preparation of the teeth.

1. Use articulating paper to **check** that there is enough room for the fibre frame and composite in occlusion. **Measure** the length

of fibre required to make a fibre frame using, for example, a periodontal probe or dental floss. Open the foil package and use tweezers to pull only the required amount of silicone bedding out of the package.

Note:

Use powder-free gloves when handling everStickC&B fibre.

2. **Cut** the required amount of fibre together with its silicone bedding, using sharp scissors. Shield the piece of fibre from light by placing it under a cover during preparation of the teeth to be bonded. Close the foil bag tightly with its sticker. Store the bag in a refrigerator (2-8°C, 35,6-46,4°F) when you do not need it.
3. **Clean** the unprepared areas of the teeth to be bonded using pumice and water, rinse with water, and air-dry the area.
4. **Etch** the surfaces of the teeth in the bonding

areas with ortho-phosphoric acid according to the bonding agent manufacturer's instructions. The recommended enamel etching time for surface-retained areas is 45 to 60 seconds. The area to be etched must be sufficiently large. Preferably etch a slightly wider area than needed rather than too small. Careful etching and bonding techniques ensure a reliable bonding of the bridge to the teeth. Rinse with water and air-dry the tooth surfaces carefully after etching.

Note: Keep the working area dry during bonding and layering of composite. Use of a rubber dam is highly recommended to keep working conditions dry.

5. **Bond** using the composite bonding technique described in the manufacturer's instructions of the bonding agent. Apply the bonding agent to the entire area to be

bonded. Light-cure the bonding agent as recommended by the manufacturer.

6. **Apply a thin layer of flowable composite** (for example, G-ænial Universal Flo) to the bonded surfaces of the teeth and/or the floor of the cavity. Do not light-cure the composite during this phase.
7. **Place the fibre.** Remove the white cover paper and use tweezers to pick the fibre up from the silicone groove. Place the fiber bundle on the tooth on top of the uncured, flowable composite. Press the other end of the fibre bundle or the entire fibre bundle tightly onto the surface of the tooth using a StickREFIX D silicone instrument, a StickSTEPPER hand instrument, or into the cavity using a StickCARRIER hand instrument. Light-cure for 5 to 10 seconds, but at the same time protect the other end of the fibre bundle

from pre-mature light-curing using a wide StickSTEPPER instrument. Position and light-cure the rest of the fibre one tooth/pontic at a time, as above. Only light-cure the fibre for 5 to 10 seconds per tooth at this stage.

Note:

- a) Spread the fibre bundle wider on the surface-retained areas to create a more extensive bonding area.
- b) Do not place the fibre too close to the gingiva, so that the cleaning spaces are not covered.
- c) Bend the fibre at the pontic area to a curved form as close as possible to the gingiva, so that maximum reinforcement is achieved. However, leave approximately 1 to 2 mm of space between the fibre and the gingiva for the composite.
- d) Do not place the fibre too close to the

gingival approxinally, so that the cleaning spaces are not covered.

- e) Position and light-cure any additional fibres as described above. Use flowable composite to attach the fibres together.
- f) When necessary, add transverse fibres to support the composite pontic.
8. **Cover and Light-Cure the fibers.** After positioning and light-curing, cover the entire fibre frame with a thin layer of composite. Light-cure the entire fibre frame for 40 seconds, one unit at a time.

Note: The fibres must be covered entirely by composite, including the interproximal areas. However, there should be enough space to allow the patient to clean the bridge and the approximal areas.
9. **Layer the pontic** using composite according to the composite manufacturer's instructions. If you do not use a rubber

dam, you can use, for example, a rubber dam strip, a sectional matrix, or a plastic strip as a moisture barrier against the gingiva. Build the core and base using flowable composite (for example, G-ænial Universal Flo). This simplifies formation of the base of the pontic. The pontic region next to the gingiva should have a light point contact, and the form should be self-cleaning. Build the dentine parts of the tooth with dentine shades and the enamel parts with enamel shades. If you wish, you can also use characterisation shades.

Note:

- a) At the location of the occlusal contact, the optimal thickness of composite to be layered on top of the fibre frame is 1-2 mm, so that it will not fracture from the fibre.

- b) Remember to preserve the cleaning spaces.

10. **Finish the bridge and adjust it into occlusion**

Note:

- a) If, after placing the fibre, you notice that the fibre is too long, shorten it using a diamond bur during the finishing phase of the bridge. Apply enamel bonding agent (for example, StickRESIN) to the exposed fibre surface to activate it and carefully remove excess bonding agent with air. Light-cure the bonding agent for 10 seconds and carefully cover the fibre again with composite.
- b) Be careful not to cut the fibres in the finishing phase

TIPS AND RECOMMENDATIONS

- The fibres must always be covered entirely with composite.

- Use a rubber dam to keep the working area dry.
- Use powder-free gloves when handling the fibres.
- During initial curing, protect the other end of the fibre to prevent premature light-curing using a StickSTEPPER instrument.
- At the occlusal contact, the optimal thickness of the composite layer on top of the fibre frame is 1-2 mm.

III CEMENTING FIBRE CONSTRUCTIONS MADE IN A DENTAL LABORATORY OR AT THE CHAIRSIDE ON A MODEL

Preparing the prosthetic work

1. Check that the fibres are visible on the cementing surfaces.

Note: The fibres must be visible on the

cementing surfaces of the work, so that the unique interpenetrating polymer network (IPN) feature of the fibres is utilised to create reliable bonding. This is especially important on surface-retained areas.

2. Remove any temporary restorations and check that the work fits.
3. Prosthetic pre-treatment
 - a) Use a carborundum bur to lightly roughen the surfaces to be cemented. Rinse with water and air-dry the surfaces.

Note: Do not use sandblasting with StickNET fibre.
 - b) Apply enamel bonding agent (for example, StickRESIN) to the roughened bonding surfaces** to activate it, protect it from light and let it take effect for 3 to 5 minutes (you can use, for instance, a metal cup as protection against light).

Carefully remove excess bonding agent with air, because an overly thick layer of bonding agent prevents the work from fitting perfectly. Light-cure the bonding agent for 10 seconds before cementation.

Note: The bonding agent used for activating the cementing surface of the fibre construction must be monomer-based and it must not contain solvents (acetone, alcohol, water). The bonding agents in the composite cement package are not necessarily suitable for activating the cementing surfaces of the fibre work

Preparing the teeth

4. Clean the surface-retained areas using a pumice and water mix.
5. Etch the surfaces of the teeth over a wide area according to the cement manufacturer's instructions. The recommended enamel etching time for surface-retained

areas is 45 to 60 seconds using a 37% orthophosphoric acid. Rinse with water and air-dry the tooth surfaces thoroughly.

6. Bond the teeth according to the cement manufacturer's instructions.

Note: Always, when possible, use a rubber dam to keep the working area dry

Cementing

7. Apply a dual cure or chemical cure composite cement onto the cementing surfaces of the work, and position the work in place.

Note: Use dual cure or chemical cure composite cements for cementing fibre work. Phosphate and glass ionomer cements are NOT suitable for cementing fibre work.
8. Remove excess cement and apply oxygen-blocking gel (for instance, glycerol gel) onto the marginal areas.
9. Light-cure dual cure cement according to the cement manufacturer's instructions.

10. Check and adjust the occlusion. Finish. Be careful not to cut the fibres when finishing the approximal areas.

IV DIRECT FIBER-REINFORCED TEMPORARY BRIDGE WITH everStickC&B

These instructions for use are intended for reinforcement of temporary bridges made with temporary crown and bridge material by using everStickC&B fibre.

1. Before preparation or extraction of the tooth/ teeth, take an impression of the working area.
2. Prepare the abutment teeth. To prepare the everStickC&B fibre frame, measure the length of fibre needed with, for example, dental floss.
3. Cut the required length of fibre together with its silicone bedding, using sharp

scissors. Shield the piece of fibre from light by placing it under a cover until it will be used. Close the foil bag tightly with its sticker, and store the bag in a refrigerator (2-8°C, 35,6-46,4°F).

4a-4d

Apply a thin layer of light curing adhesive* (e.g. StickRESIN) on top of the abutments. Do not light-cure the adhesive yet. Place the fibre bundle on top of the abutments. Use a StickREFIX D silicone instrument to press the fibre into the correct position and shape. To achieve maximum reinforcement, bend the fibre in the pontic area to a curved form as close as possible to the gingiva. However, leave approximately 1-2 mm of space between the fibre and the gingiva for the temporary crown and bridge material. Aproximally do not place the fibre too close to the gingiva, so as not to cover

the cleaning spaces. Initially light-cure the fibre through the silicone for 10 seconds. Finalise the polymerisation by light-curing the entire fibre frame for 40 seconds per area. *** NOTE:** The light curing adhesive used must be monomer-based and must not contain etchants or solvents (acetone, alcohol, water, etc.). Do not etch the abutments or use dentin primer prior to adhesive application.

TIP: If there is enough space for extra fibres, you can add more fibre layers to increase the strength of the bridge.

5. Follow the user instructions of the temporary crown and bridge material manufacturer. Dispense temporary material first between the bridge frame and the gingiva, to eliminate air gaps. Fill the impression with temporary material, place it over the prepared teeth, and allow it to set.

6. When the temporary material has set, remove the bridge, along with the impression material.
7. Contour and polish the temporary bridge according to the manufacturer's instructions and cement it with suitable temporary cement. Check and adjust the occlusion.

STORING: everStick products should always be stored in a refrigerator (2-8°C, 35,6-46,4°F). In addition, the products should be protected from light by packing them in the sealed foil package after use. An elevated temperature and exposure to bright light may shorten the lifetime of everStick products.

Prior to application, the products are taken out of the refrigerator and the foil package opened, but kept away from bright daylight or artificial light. While cutting the fibre bundle, the rest of the fibre bundle inside the foil package should

be kept covered from light. Immediately after cutting a sufficient length for the fibre construction, the foil package is carefully resealed and returned to the refrigerator. (Shelf life: 2 years from date of manufacture)

PACKAGES

everStick Intro

8 cm everStickC&B; 8 cm everStickPERIO;
1x StickREFIX D silicone instrument

everStick COMBI

8cm everStickC&B; 8cm everStickPERIO;
30 cm² everStickNET; 5x everStickPOST Ø 1.2;
1x StickSTEPPER hand instrument;
1x StickREFIX D silicone instrument

everStick Starter Kit:

8 cm everStickC&B; 8 cm everStickPERIO;

30 cm² everStickNET; 5 x everStickPOST
Ø 1.2; 5 ml StickRESIN bottle; 2 ml G-ænial
Universal Flo syringe A2; 20 dispensing tips,
1 light protective cover; 1x StickREFIX D
silicone instrument; 1x StickSTEPPER;
1x StickCARRIER

Refills

2 x 12 cm fibre bundle
1 x 8 cm fibre bundle

NOTE: The everStick products should be used clinically with care and the patient should be warned not to abrade the fitting surface so as to avoid exposing irritation-causing fibres. If the surface of the fibre bundle feels dry, but it is fully bendable and not polymerized, adding a drop of resin (such as StickRESIN) will return the flexibility/workability of the material.

Polymerization can be observed as white spots at bending area when bending the bundle.

The everStick fibres do not achieve their full strength immediately after the final light-curing of 40 seconds. The polymerization of the fibres will still continue during the next 24 hours.

StickSTEPPER, StickCARRIER hand instruments and StickREFIX D, StickREFIX L silicone instruments must be sterilized before use.

WARNING: Personal protective equipment (PPE) such as gloves, face masks and safety eyewear should always be worn. Unpolymerised resin can cause skin sensitisation to acrylates in some people. If your skin comes in contact with resin, wash it thoroughly with soap and water. Avoid contact of uncured material with skin, mucous membrane, or eyes. Unpoly-

merised everStick products may have a slight irritating effect and lead to sensitization to methacrylates in rare cases. The use of powder free gloves is recommended with everStick products. Polymerize everStick before waste disposal.

US Federal law restricts this device to sale by or on the order of a dentist.

Some products referenced in the present IFU may be classified as hazardous according to GHS. Always familiarize yourself with the safety data sheets available at:
<http://www.gceurope.com>
or for The Americas
<http://www.gcamerica.com>
They can also be obtained from your supplier.

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