


RELIABILITY WITH PATENTED ACTIVE-CONNECT-TECHNOLOGY (ACT)

STANDARD DUAL-CURING LUTING COMPOSITES

Dual-curing composites without ACT often do not cure completely on acid-etched tooth surfaces and/or those pretreated with adhesive or primer and require an additional activator (referred to as self-cure activators (SCA)).

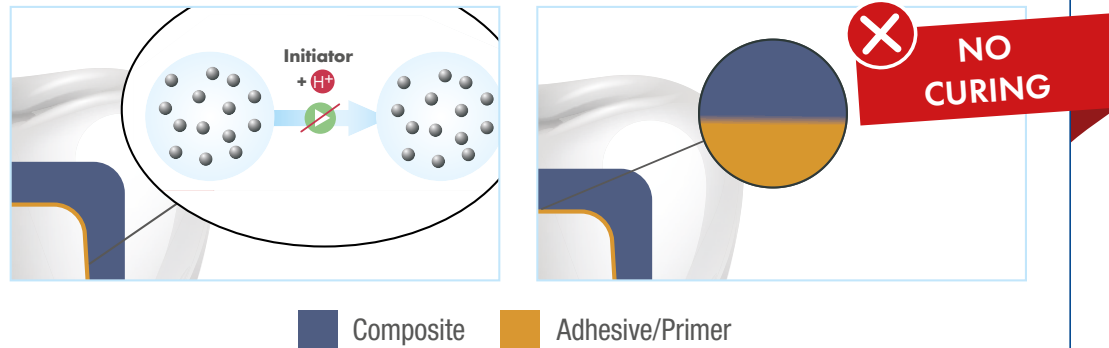

Please contact us if you have any more questions about the adhesive luting composite Visalys® CemCore:

Phone: **+49 (0) 27 74 / 705 0** or


E-mail: **visalys-cemcore@kettenbach.com**

Kettenbach GmbH & Co. KG · Im Heerfeld 7 · 35713 Eschenburg · www.kettenbach-dental.com

VISALYS® CEMCORE WITH ACT



The integrated phase transfer catalysts lead to mixing at the interface between Visalys® CemCore and Visalys® Tooth Primer. The integrated AC Activator guarantees reliable polymerization, even under acidic conditions or in the dark.

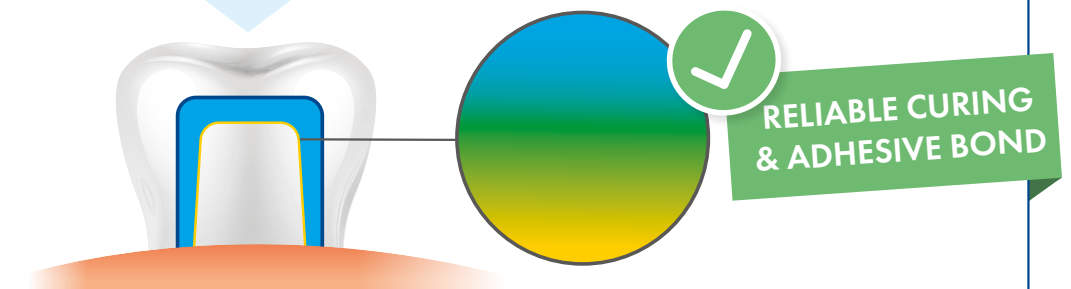

Composite (Visalys® CemCore) Primer (Visalys® Tooth Primer)

Dieses Dokument wird Ihnen zur Verfügung gestellt von:

Med-Dent24 Handelsgesellschaft mbH
Gotenweg 7
63128 Dietzenbach

Telefon: +49 (0)6074 803 3615
Telefax: +49 (0)6074 803 4092
E-Mail: info@med-dent24.com

Registergericht: Offenbach am Main
Registernummer: HRB 56518
Umsatzsteuer-Identifikation-No.: DE307869540

Geschäftsführerin: Ana Ban